Instruction for use

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**A SOLID-PHASE ENZYME IMMUNOASSAY
FOR THE QUANTITATIVE DETERMINATION
OF ALPHA-FETOPROTEIN
IN HUMAN SERUM OR PLASMA**

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
Authorized Representative in EU:
Polmed.de
Steinacker 20, D-73773
Aichwald, Germany
e-mail: info@polmed.de

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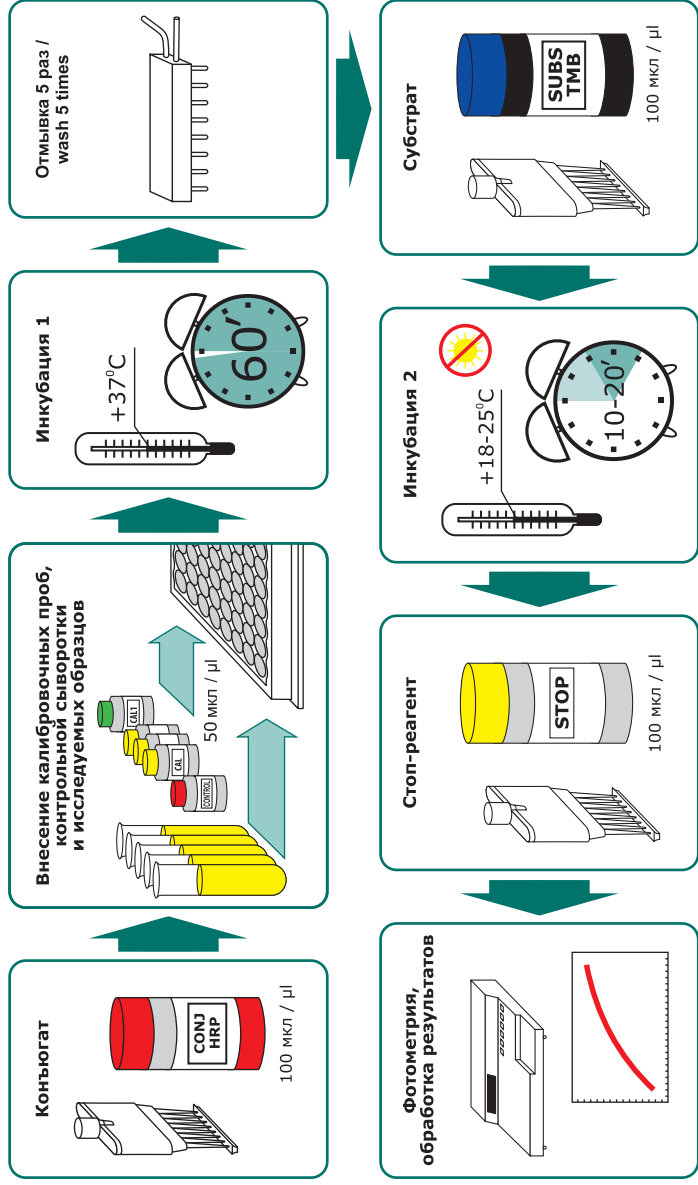CONTENT

1. INTENDED USE	10
2. SUMMARY AND EXPLANATION	10
3. PRINCIPLE OF THE TEST	10
4. WARNINGS AND PRECAUTIONS	11
5. KIT COMPONENTS	12
6. SPECIMEN COLLECTION AND STORAGE	13
7. TEST PROCEDURE	13
8. QUALITY CONTROL	15
9. CALCULATION OF RESULTS	15
10. EXPECTED VALUES	15
11. PERFORMANCE CHARACTERISTICS	16
12. LITERATURE	16

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1. Johnson, P.J, (2002) Tumor Markers in Primary Malignancies of the liver. In «Tumor Markers: Physiology, pathobiology, technology and clinical applications», ed. Dimandis E.P. AACCC Press, Washington pp 269-276.
2. Stenman, U-H and Alfthan, H. (2002) Markers for Testicular Cancer. In «Tumor Markers: Physiology, pathobiology, technology and clinical applications», ed. Dimandis E.P. AACCC Press, Washington, pp 351-359.
3. Christiansen, M. et al. Alpha-fetoprotein in plasma and serum of healthy adults: preanalytical, analytical and biological sources of variation and construction of age-dependent reference intervals. Scand J Invest 2001 61: 205-216.
4. Trape, J. et al. Reference change value for a-Fetoprotein and its application in early detection of hepatocellular carcinoma in patients with hepatic disease. Clin Chem 2003 49(7): 1209-1211.

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Instruction for use

A SOLID-PHASE ENZYME IMMUNOASSAY FOR THE QUANTITATIVE DETERMINATION OF ALPHA-FETOPROTEIN IN HUMAN SERUM OR PLASMA

1. INTENDED USE

A solid-phase enzyme immunoassay for the quantitative determination of alpha-fetoprotein in blood serum or plasma.

This kit is designed for measurement of alpha-fetoprotein in blood serum or plasma. For possibility of use with other sample types, please, refer to Application Notes (on request). The kit contains reagents sufficient for 96 determinations and allows to analyze 42 unknown samples in duplicates.

2. SUMMARY AND EXPLANATION

Alpha-fetoprotein (AFP) is a glycoprotein with a MW ca. 65 kDa which is secreted by fetal liver and yolk sac. AFP represents the main protein of fetal serum while being found in trace quantities in adults. Serum AFP quantitative determination is used in primary diagnostics and monitoring of hepatocellular liver cancer, trophoblastic tumours of testicles and ovary as well as theratomas and theratocarcinomas. Quantitative determination of AFP in serum of pregnant women or in amniotic fluid during week 15-20 of gestation is widely used for laboratory screening of Down syndrome and defects of spinal cord.

3. PRINCIPLE OF THE TEST

This test is based on two-site sandwich enzyme immunoassay principle. Tested specimen is placed into the microwells coated by specific murine monoclonal to human AFP-antibodies. Antigen from the specimen is captured by the antibodies coated onto the microwell surface. Second antibodies – murine monoclonal to human AFP, labelled with peroxidase enzyme, are then added into the microwells. After washing procedure, the remaining enzymatic activity bound to the microwell surface is detected and quantified by addition of chromogen-substrate mixture, stop solution and photometry at 450 nm. Optical density in the microwell is directly related to the quantity of the measured analyte in the specimen.

4. WARNINGS AND PRECAUTIONS

4.1. For professional use only.

4.2. This kit is intended for in vitro diagnostic use only.

4.3. INFECTION HAZARD: There is no available test methods that can absolutely assure that Hepatitis B and C viruses, HIV-1/2, or other infectious agents are not present in the reagents of this kit. All human products, including patient samples, should be considered potentially infectious. Handling and disposal should be in accordance with the procedures defined by an appropriate national biohazard safety guidelines or regulations.

4.4. Avoid contact with stop solution containing 5.0% H₂SO₄. It may cause skin irritation and burns.

4.5. Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents may give false results.

4.6. Do not use the kit beyond the expiration date.

4.7. All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microplate readers.

4.8. Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.

4.9. Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guidelines or regulations.

4.10. Do not mix reagents from different lots.

4.11. Replace caps on reagents immediately. Do not swap caps.

4.12. Do not pipette reagents by mouth.

4.13. Specimens must not contain any AZIDE compounds – they inhibit activity of peroxidase.

4.14. Material Safety Data Sheet for this product is available upon request directly from XEMA Co., Ltd.

4.15. The Material Safety Data Sheet fit the requirements of EU Guideline 91/155 EC.

5.1. Contents of the Kit

5. KIT COMPONENTS

Symbol	Description	Qty	Units	Colour code	Stability of opened/diluted components	
1	SORB MTP AFP EIA strips, 8x12 wells	polystyrene microwells coated with murine monoclonal to human AFP	1	pcs	until exp.date	
2	CAL 1-6 Calibrator set, 0.8 ml each, zero calibrator C1 - 6 ml. The set contains 6 calibrators: 0; 5 ;15; 50; 150; 500 IU/ml	human alpha-fetoprotein diluted in tris buffered BSA solution, preservative – 0.01% Bronidox L, 0.01% 2-Methyl-4-isothiazolin-3-one-hydrochloride; also contains red dye	6	pcs	red (C1 – colourless)	2 months
3	CONTROL Control serum (0.8 ml)	dilution of preselected human serum, with high content of alpha-fetoprotein with BSA solution; preservative – 0.01% Bronidox L, 0.01% 2-Methyl-4-isothiazolin-3-one-hydrochloride, colourless	1	pcs	colourless	2 months
4	CONJ HRP Conjugate, 14 ml	aqueous solution of murine monoclonal to human AFP coupled with horseradish peroxidase diluted on phosphate buffered solution with casein from bovine milk and detergent (Tween-20), contains 0.1% phenol as preservative and red dye	1	pcs	red	until exp.date
5	SUBS TMB Substrate solution, 14 ml	ready-to-use single-component tetramethylbenzidine (TMB) solution	1	pcs	colourless	until exp.date
6	BUF WASH 26X Washing solution concentrate 26X, 22 ml	aqueous solution of sodium chloride and detergent (Tween 20), contains proClin300 as a preservative	1	pcs	colourless	Concentrate – until exp.date Diluted washing solution – 45 days at 2-8 °C or 15 days at RT
7	STOP Stop solution, 14 ml	5.0% vol/vol solution of sulphuric acid	1	pcs	colourless	until exp.date
8	N003 Plate sealing tape		2	pcs		N/A
9	K2251 Instruction AFP EIA		1	pcs		N/A
10	K225Q QC data sheet AFP EIA		1	pcs		N/A

5.2. Equipment and material required but not provided

- Distilled or deionized water;
- Automatic or semiautomatic multichannel micropipettes, 100–250 µl, is useful but not essential;
- Calibrated micropipettes with variable volume, range volume 25–250 µl;
- Dry thermostat for 37 °C ±0.1 °C;
- Calibrated microplate photometer with 450 nm wavelength and OD measuring range 0–3.0.

5.3. Storage and stability of the Kit

Store the whole kit at +2...+8 °C upon receipt until the expiration date.

After opening the pouch keep unused microtiter wells TIGHTLY SEALED BY ADHESIVE TAPE (INCLUDED) to minimize exposure to moisture.

6. SPECIMEN COLLECTION AND STORAGE

This kit is intended for use with serum or plasma (ACD- or heparinized). Grossly hemolytic, lipemic, or turbid samples should be avoided.

Specimens may be stored for up to 48 hours at +2...+8 °C before testing. For a longer storage, the specimens should be frozen at -20 °C or lower. Repeated freezing/thawing should be avoided.

7. TEST PROCEDURE**7.1. Reagent Preparation**

- All reagents (including unsealed microstrips) should be allowed to reach room temperature (+18...+25 °C) before use.
- All reagents should be mixed by gentle inversion or vortexing prior to use. Avoid foam formation.
- It is recommended to spin down shortly the tubes with calibrators on low speed centrifuge.
- Prepare washing solution from the concentrate BUF WASH 26X by 26 dilutions in distilled water.

7.2. Procedural Note:

It is recommended that pipetting of all calibrators and samples should be completed within 3 minutes.

7.3. Assay flowchart

See the example of calibration graphic in Quality Control data sheet.

7.4. Alternative units:

1 IU/ml = 1.25 ng/ml

7.5. Assay procedure

1	Put the desired number of microstrips into the frame; allocate 14 wells for the calibrators CAL 1-6 and control samples CONTROL and two wells for each unknown sample. DO NOT REMOVE ADHESIVE SEALING TAPE FROM UNUSED STRIPS.
2	If suggested analyte concentration in the sample exceeds the highest calibrator, additionally dilute this sample accordingly, using zero calibrator. Use of other buffers or reagents for sample dilution may lead to incorrect measurement.
3	Dispense 100 µl of CONJ HRP into the wells.
4	Pipet 50 µl of calibrators CAL 1-6, control samples CONTROL and unknown samples into the wells. Cover the wells by plate adhesive tape (included into the kit).
5	Incubate 60 minutes at 37 °C.
6	Prepare washing solution by 26X dilution of washing solution concentrate (BUF WASH 26X) by distilled water. Wash the strips 5 times.
7	Dispense 100 µl of SUBS TMB into the wells.
8	Incubate 10-20 minutes at 18...+25 °C.
9	Dispense 100 µl of STOP into the wells.
10	Measure OD (optical density) at 450 nm.
11	Set photometer blank on first calibrator.
12	Apply point-by-point method for data reduction.

7.6. Handling notes

Calibrators and control sample(s) – only one freezing/thawing cycle is allowed.

8. QUALITY CONTROL

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results.

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable federal, state, and local standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.

The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications.

9. CALCULATION OF RESULTS


9.1. Calculate the mean absorbance values (OD450) for each pair of calibrators and samples.

9.2. Plot a calibration curve on graph paper: OD versus alpha-fetoprotein concentration.

9.3. Determine the corresponding concentration of alpha-fetoprotein in unknown samples from the calibration curve. Manual or computerized data reduction is applicable on this stage. Point-by-point or linear data reduction is recommended due to non-linear shape of curve.

9.4. Below is presented a typical example of a standard curve with the XEMA Co. Not for calculations!

Calibrators	Value	Absorbance Units (450 nm)
CAL 1	0 IU/ml	0.05
CAL 2	5 IU/ml	0.098
CAL 3	15 IU/ml	0.604
CAL 4	50 IU/ml	1.334
CAL 5	150 IU/ml	1.893
CAL 6	500 IU/ml	2.230

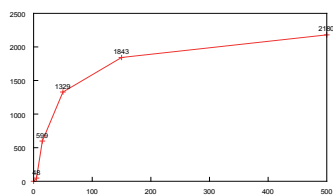10. EXPECTED VALUES

Therapeutical consequences should not be based on results of IVD methods alone – all available clinical and laboratory findings should be used by a physician to elaborate therapeutically measures. Each laboratory should establish its own normal range for AFP. Based on data obtained by XEMA, the following normal range is recommended (see below). NOTE: the patients that have received murine monoclonal antibodies for radioimaging or immunotherapy develop high titered anti-mouse antibodies (HAMA). The presence of these antibodies may cause false results in the present assay. Sera from HAMA positive patients should be treated with depleting adsorbents before assaying.

Sex, age	Units, IU/ml		Units alternative, ng/ml	
	Lower limit	Upper limit	Lower limit	Upper limit
Healthy donors	-	10.0	-	12.5

Pregnancy, week	Median, IU/ml	SKO
14	21.7	0.43
15	28.3	0.47
16	30.0	0.51
17	36.3	0.49
18	43.8	0.52
19	53.3	0.50
20	60.0	0.55
21	63.3	0.57

11. PERFORMANCE CHARACTERISTICS

11.1. Analytical specificity / Cross reactivity

Analyte	Cross-reactivity, % wt/wt
Albumin	<0.1
hCG	<0.1
Lactogen	<0.1

11.2. Analytical sensitivity

Sensitivity of the assay was assessed as being 0.9 IU/ml.

11.3. Linearity

Linearity was checked by assaying dilution series of 5 samples with different alpha-fetoprotein concentrations. Linearity percentages obtained ranged within 90 to 110%.

11.4. Recovery

Recovery was estimated by assaying 5 mixed samples with known alpha-fetoprotein concentrations. The recovery percentages ranged from 90 to 110%.

11.5. Hook-effect

No hook-effect has been noticed with samples up to 12000 IU/ml.

12. LITERATURE

1. Johnson, P.J, (2002) Tumor Markers in Primary Malignancies of the liver. In «Tumor Markers: Physiology, pathobiology, technology and clinical applications», ed. Dimandis E.P. AACC Press, Washington pp 269-276.
2. Stenman, U-H and Alfthan, H. (2002) Markers for Testicular Cancer. In «Tumor Markers: Physiology, pathobiology, technology and clinical applications», ed. Dimandis E.P. AACC Press, Washington, pp 351-359.
3. Christiansen, M. et al. Alpha-fetoprotein in plasma and serum of healthy adults: preanalytical, analytical and biological sources of variation and construction of age-dependent reference intervals. Scand J Invest 2001 61: 205-216.
4. Trape, J. et al. Reference change value for a-Fetoprotein and its application in early detection of hepatocellular carcinoma in patients with hepatic disease. Clin Chem 2003 49(7): 1209-1211.

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Russian Diagnostic
Manufacturers Association

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Russian Association
of Medical Laboratory
Diagnosticians

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